

# IMAGEnet i-base

## User manual

This manual is based on IMAGEnet i-base advanced.  
It will be mentioned if other licenses are required.

**Topcon Europe Medical BV**

Version 1.16

2017-01-09

 **0344**

## CONTENTS

<b>CONTENTS .....</b>	<b>1</b>
<b>PREFACE .....</b>	<b>2</b>
INTENDED USE .....	2
DISPLAYS FOR SAFE USE .....	2
SAFETY CAUTIONS .....	2
USER MAINTENANCE .....	3
MEASUREMENTS .....	4
<b>START IMAGENET I-BASE .....</b>	<b>6</b>
LOGIN.....	6
AGREE LICENSE AGREEMENT.....	7
RISK CONTROLS RELATED TO THE INSTRUCTIONS FOR USE .....	8
<b>IMAGENET I-BASE LAYOUT'S .....</b>	<b>9</b>
CREATE LAYOUT.....	9
SELECT A LAYOUT.....	10
CONFIGURE A STARTUP LAYOUT .....	10
SHORTCUTS WINDOW .....	10
<b>CAPTURE SESSION PROCEDURE (ONLY WITH EYEDOCTOR LICENSE) .....</b>	<b>13</b>
SELECT CAPTURE LAYOUT.....	13
START A CAPTURE SESSION.....	13
ADD A NEW PATIENT OR SELECT A PATIENT.....	13
ADD A NEW PATIENT.....	14
SELECT A PATIENT .....	15
CAPTURE WINDOW READY TO CAPTURE IMAGES.....	16
START TIMER .....	16
DELETE IMAGES DURING CAPTURE .....	17
UNDELETE IMAGES DURING CAPTURE .....	17
SAVE CAPTURED IMAGES .....	18
<b>PATIENT MANAGEMENT / DATA .....</b>	<b>19</b>
ADD NEW PATIENT .....	19
EDIT (MANAGE) PATIENT INFORMATION .....	20
SHOW PATIENT DATA .....	21
OPEN PATIENT DATA.....	22
WORKING IN A MULTI-UP .....	23
COMPOSITION TOOL BAR .....	24
FITTING IMAGES BY ZOOM .....	25
IMAGE PROCESS .....	26
TOOLS WINDOW .....	27
PLUGIN WINDOW .....	28
MODULES .....	29
MEASUREMENTS .....	33
PRINT PATIENT DATA .....	34
PRINT PREVIEW .....	35

## PREFACE

This manual is valid for IMAGEnet i-base version 3.18.4 and higher. For more information see the installation manual. Electronic instructions for use are published in PDF format. A PDF reader can be downloaded at [www.adobe.com](http://www.adobe.com). This instructions for use can be requested on paper if required.



## INTENDED USE

With IMAGEnet i-base, data of instruments can be acquired, analysed, printed and stored in the database along with relevant patient information. Supported data includes images, movies and refraction data. IMAGEnet i-base allows the user to do measurements and enhancements on ophthalmology images.

## DISPLAYS FOR SAFE USE

To encourage safe and proper use and to prevent loss of data, important messages are inserted in the User Manual.

Presentation of notes and cautions in this manual are done confirm the following table:

Icon		Meaning
	<b>NOTE</b>	This symbol indicates that the message has special interest or importance to the reader. Failure to read the note will not result in physical harm to the reader, equipment or data.
	<b>CAUTION</b>	This symbol indicates that if this information is ignored, data loss or material damage may occur.

## SAFETY CAUTIONS



### CAUTION

Please read the user manual carefully before using IMAGEnet i-base.



### NOTE

We strive for a manual with screenshots of the latest version of IMAGEnet i-base, but it could be possible that there are screenshots of older versions.



### NOTE

This user manual is written on the assumption that the user already has knowledge of the personal computer (PC) and Microsoft Windows operating systems. This user manual does not describe the use of a PC or Microsoft Windows.



### NOTE

Should you have unexpected problems, please contact your local Topcon dealer.

**CAUTION**

Do not use IMAGEnet i-base when there is visible damage to cables connecting the PC to the instrument or visible damage to any other part of your set up. In such event please contact your local Topcon representative.

**CAUTION**

Make sure that you always have an up to date virus scanner and firewall installed on your system.

**CAUTION**

Avoid disconnecting cables from the PC when the IMAGEnet i-base software is running.

**CAUTION**

Avoid resetting or turning off the PC when the IMAGEnet i-base software is running.

**CAUTION**

Avoid resetting or turning off the instrument when the IMAGEnet i-base software is running.

**CAUTION**

A diagnosis cannot be based only on measurements or methods from IMAGEnet i-base. The doctor should consult multiple sources before giving a diagnosis.

**CAUTION**

Do not use IMAGEnet i-base when there are indications of problems or damage to the system. In such event please contact your local Topcon representative.

## USER MAINTENANCE

### UNEXPECTED EVENTS

In case of an unexpected event, such as a power failure, PC crash or i-base crash, IMAGEnet i-base is shut down without saving, the following actions need to be taken.

When possible, restart i-base in the usual way. i-base will open all data that was buffered at the time of the event and display the information on the data that was available in the capture tab. Additional information will be available from the data properties window.

This data needs to be reviewed and stored at the correct patient.

If, due to circumstances, the patient to which the data belongs cannot be deduced, this data must absolutely not be used.

Should you be unable to do this, or have different problems, please contact your local Topcon representative.

### STORAGE LOCATION

When running IMAGEnet i-base as a capture station in a network, it is recommended to create an extra database and image save location on the local PC where i-base is installed. Should you have any problems with your network connection, you will be able to switch to this location to save your data. This data needs to be moved to your normal database when your network connection becomes available again.

To do so create an extra database in the usual way. In case of problem you can switch to this database by logging in (select File > Login). In the dialog box: select your local database and enter your credentials and click Ok. Also make sure that your local save location is selected (File > Options > File locations > Current image save location). Your data will now be saved locally.

When you wish to save in the usual location again repeat these steps for your network save location. Data stored locally needs to be moved to the usual location using Patient packs, see the Help file: Patient packs.

## MEASUREMENTS



### CAUTION

For maximal accuracy in measurements made on fundus images the Littman formula needs to be used to correct the measured values for the aberrations in the human eye.

## MEASUREMENTS IN FUNDUS IMAGES

The accuracies for all measurements are estimated to be 10 %



### CAUTION

The accuracy of line or area measurements conform to those of the Gullstrand test eye.

## MEASUREMENTS IN OCT DATASETS

The measurement function on the OCT images – caliper function – is enabled for Posterior Segment images and can be used for Cornea Thickness measurements on Anterior segment images.

The following table gives an overview of the applied measurement function:

	Vertical measurement	Horizontal measurement
Anterior	Applied for Cornea Thickness	Applied for Cornea Thickness
Posterior (Hand-held lens)	Applied	Not applied
Posterior (FV-1L)	Applied	Applied

### Posterior measurements

In case of Posterior Segment images the vertical axis i.e. depth axis or z axis is always enabled and independent of the optics used. The horizontal component of the measurement function strongly depends on the magnification of the optics, and on the alignment of the components in the optical system. Since the alignment can not be guaranteed for hand-held lenses, the horizontal component is only enabled when the Fundus Viewer FV-1L is used.

A combination of vertical and horizontal component allows the user to measure free-forms on the OCT images.

The system tolerances of the measurement components are:

**Vertical Measurement is a summation of pixels where the size of each individual pixel is 4.2  $\mu\text{m}$  with an error < 2%.**

Horizontal Measurement is a summation of pixels where the size of each individual pixel depends on the optical magnification. The error of the measurement component is <10% for the central region of the retina. The central region of the retina is defined as an area with a radius of approximately 4 mm around the fovea.

**CAUTION**

The horizontal component of the measurement function depends on the relation between the optics settings i.e. magnification of the slitlamp and hand-held lens AND the settings selected in the application plug-in. The horizontal measurement function is not reliable if this is not set correctly by the user.

**CAUTION**

The accuracy of the horizontal component of the measurement function is based on the Gullstrand Eye conditions.

*Anterior Top of Cornea measurement*

In case of Anterior Segment images the vertical axis i.e. depth axis or z axis is always enabled and independent for the optics used. The horizontal component of the measurement function strongly depends on the magnification of the optics, and on the alignment of the components in the optical system. Since the alignment cannot be guaranteed for all the possible variation of the position of the Slitlamp in relation to the position of the anterior segment of the eye, the horizontal component is only intended to be used as part of the measurement of the Cornea Thickness, where the error contribution of these effects are limited.

In order to perform this measurement take an image by aiming at the centre of the pupil and adjust the SL SCAN-1 alignment until you find the optical reflex, that will be clearly visible as a bright vertical line in the SL SCAN-1 image. Usually this is done by moving the slitlamp up or down a little.

When you find the optical reflex, press the trigger to take an image. You will need to save this image to be able to measure on it. Measurement will need to be done just next to the optical reflex.

Top of Cornea measurement is a summation of pixels where the size of each individual pixel depends on the optical magnification. **The error of the measurement component is < 5% for the central region of the cornea.** The central region of the cornea is defined as an area with a diameter of approximately 4 mm around the centre of the pupil.

**CAUTION**

Anterior measurement is intended for central cornea thickness measurement only.

**CAUTION**

The horizontal component of the measurement function depends on the relation between the optics settings i.e. magnification of the slitlamp AND the settings selected in the application plug-in. The measurement function is not reliable if this is not set correctly by the user.

**CAUTION**

The accuracy of the Top of Cornea measurement function is based on comparison to OCT 2000 FA plus and measurement of glass plates of known thickness.

## MEASUREMENTS IN OCT DATASETS


The measurement function on the OCT images – caliper function – is enabled for Posterior Segment images and can be used for Cornea Thickness measurements on Anterior segment images.

## START IMAGENET I-BASE

### LOGIN



To startup IMAGENet ibase double click the shortcut icon on your desktop.  
The IMAGENet ibase login window will appear.

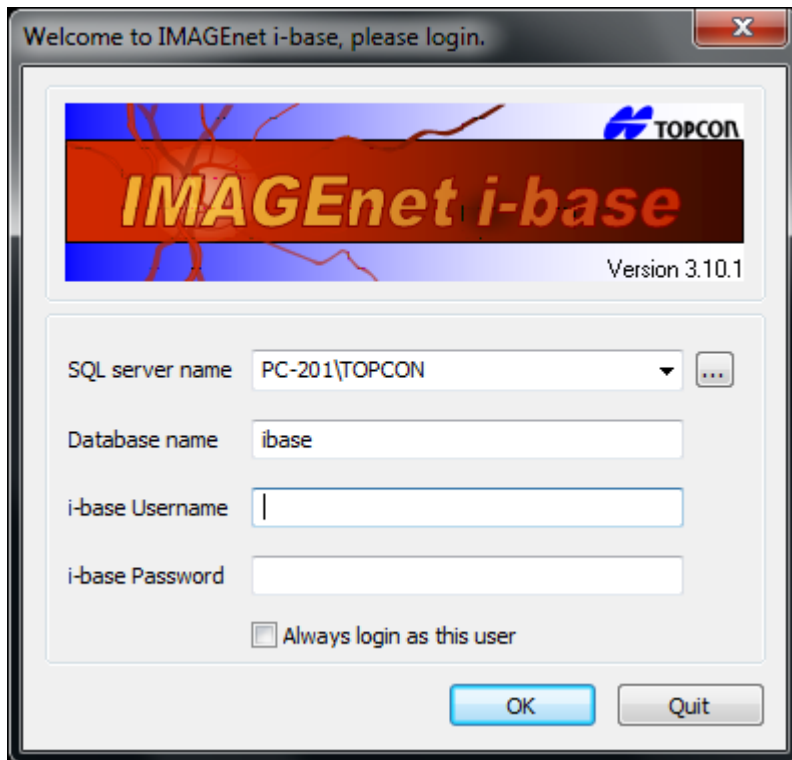

If only one IMAGENet i-base user is active on this pc, you can mark "Always login as this user"

After filling in all necessary fields click "OK"


As Default the database name is "ibase" username and password are "ibaseAdmin". If these parameters are not correct to login please contact your dealer.



#### NOTE

For safety reasons it is recommended to change this default password.  
To do this goto IMAGENet i-base menu: **File > Change password.**

## AGREE LICENSE AGREEMENT

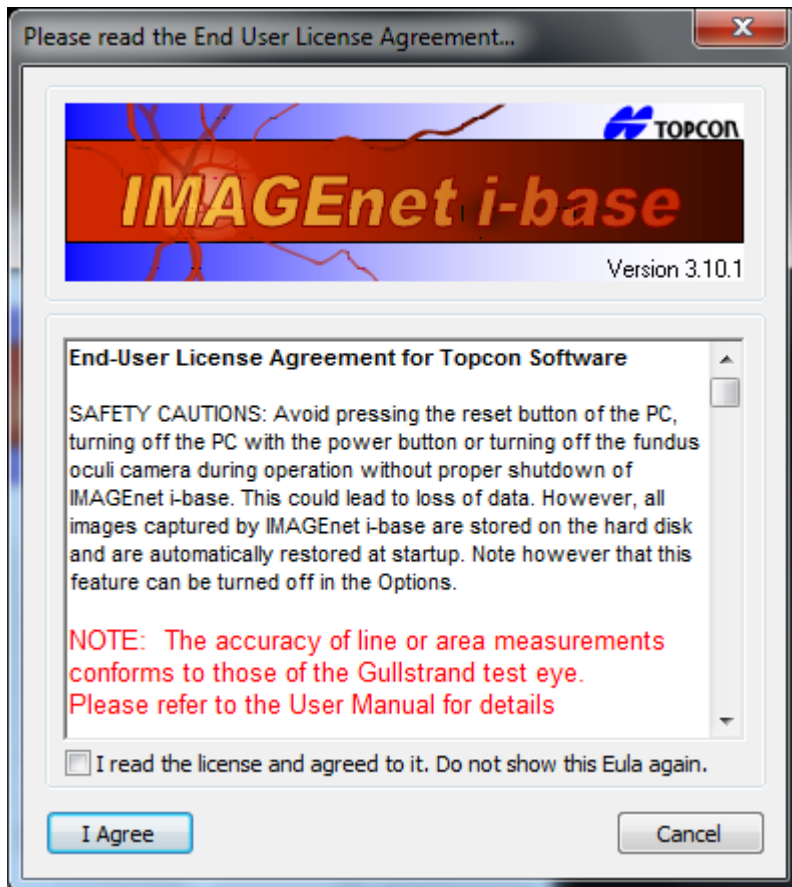


Read the end user license agreement for Topcon software and Press **"I Agree"** to continue with ibase.



## RISK CONTROLS RELATED TO THE INSTRUCTIONS FOR USE

The end-user is notified at start-up of the software with the 'End-User License Agreement' as shown below. This agreement will be shown to all (new) users until they agree to it.

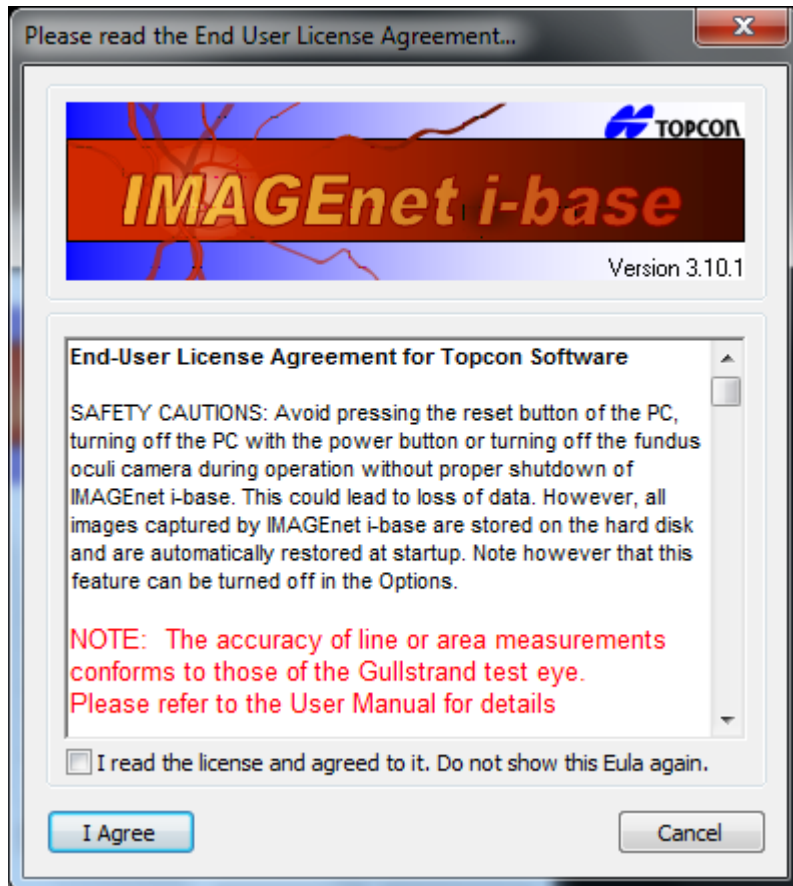


## IMAGENET I-BASE LAYOUT'S

A layout in IMAGEnet i-base is the i-base screen overview with all the shown windows and toolbars. Also i-base main window position (including which monitor ) and size (including minimized or maximized state) are parts of the layouts. Settings and options for windows are not included in a layout, except for the thumbnail window state.

From the View menu is it possible to show and hide all the windows and toolbars to create a different layout. Most bars and windows can be moved, docked and undocked with the normal Windows functionality.


It can be useful to have different layouts for different i-base actions. For example, capturing and viewing.

Layouts can be managed with the 'Manage layouts' dialog, which can be started from the menu (View > Layouts > Manage layouts) or from the 'Windows toolbar'.

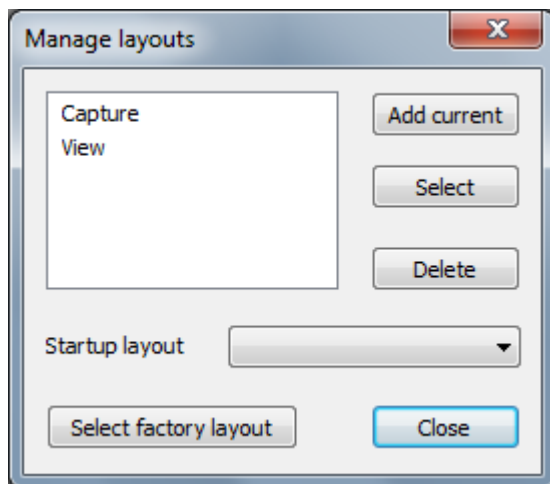

## CREATE LAYOUT

To create a new layout press in the 'Manage layouts' dialog the 'Add current' button and enter the preferred name. This saves the current layout.

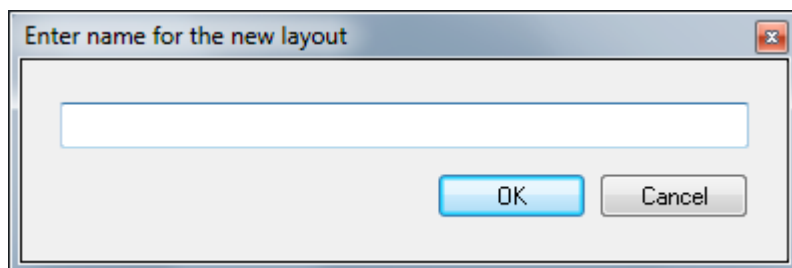

The same action can be done with the 'Save current layout' button in the 'Windows' toolbar.



## SELECT A LAYOUT

Selecting a layout can be done from the layout management by choosing the layout from the list box and pressing the 'Select' button.

Another way is to select the preferred layout from the dropdown in the 'Windows' toolbar.



## CONFIGURE A STARTUP LAYOUT

If it is desired to startup i-base always with the same layout, select in the 'Manage layout' dialog the preferred layout at 'Startup layout'.




### NOTE

If no layout is selected as 'Startup layout' IMAGEnet i-base will remember his last layout and startup with it the next time.


## SHORTCUTS WINDOW

If you want to make the buttons smaller or bigger in the Shortcuts windows, just hold down the Control key on your keyboard and click one time in an empty space on the Shortcuts window. Now scroll with your mouse wheel while still holding the control key.

To customize the shortcut bar right click on the shortcut bar and click "Customize shortcut bar".

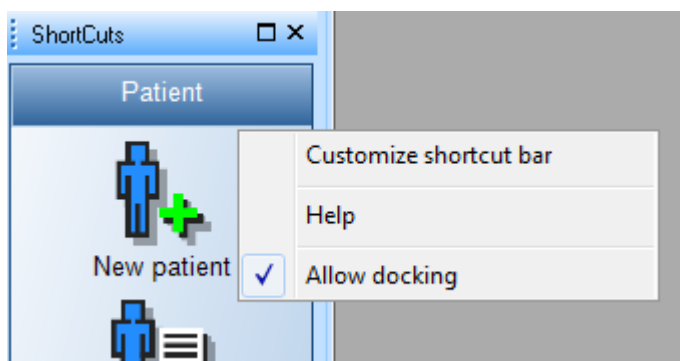

A screen appears to configure the tabs:

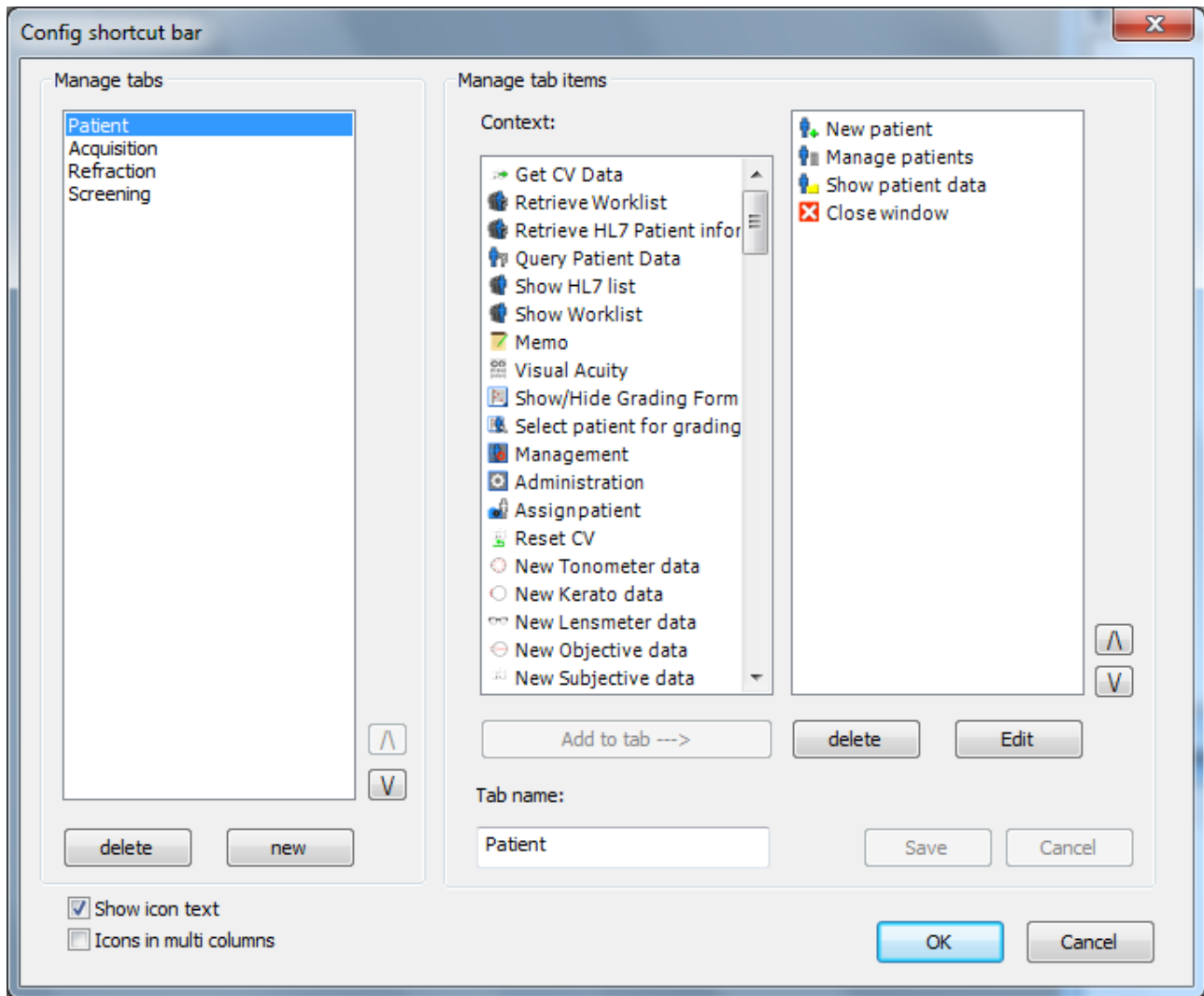


### Add a tab

To add a new tab press "new" and enter the name of the tab. Add item(s) by clicking an item and then the button "Add to tab --->" or double click the item. Press "Save" when you are satisfied with the new tab.

### Modify a tab

To modify a tab, just select the tab name you want. Now it is possible to add, delete or modify the items in this tab. Modifying of the items in tab can be done with the "Edit" button or double click. A new dialog appears where it is possible to change the text and the icon of the item. Press "Save" when you are done.

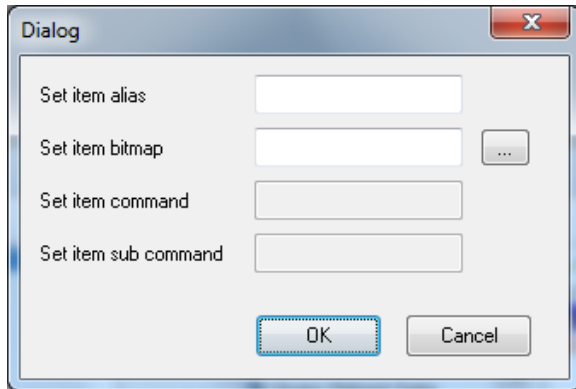

### Delete a tab/item

To delete a tab/item select the tab/item name you want to delete and press the "delete" button. The tab with all the items will be removed from the shortcut bar and the item will be removed from your tab. Press "Save" when you are done.

### Start an external file

To start from i-base an external program or open a file.

Examples:

To open a MS Word document:

alias: Open my file  
 bitmap: choose here a picture  
 commando: browse to WINWORD.EXE  
 sub commando: the file you want to open

To start Excel:

alias: Start Excel  
 bitmap: choose here a picture  
 commando: browse to EXCEL.EXE  
 sub commando: leave it empty

### Start/Config a plugin:

To start or config a plugin.

Example:

To Start the Cell Count module :

alias: Start Cell Count  
 bitmap: choose here a picture  
 commando: Cell Count  
 sub commando: Cell Count (automatically)

## CAPTURE SESSION PROCEDURE (ONLY WITH EYEDOCTOR LICENSE)


### SELECT CAPTURE LAYOUT

If there is created a special layout for capturing select this from the 'Windows' toolbar or the 'View' menu. If there is no layout for capturing show the 'Capture session control' window by selecting this window from the 'View' menu.

### START A CAPTURE SESSION

From the capture session control window you can select the different capture procedures.

In this menu you can start a new procedure by double click or select a procedure and press the "start" button with the mouse pointer.

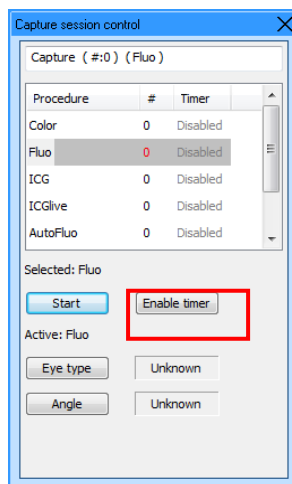

**Note,** When you select for example "Fluo" procedure, a capture session window will come up. IMAGEnet Ibase is now ready to receive images.

**Note,** To enable timer press the "Enable timer" button. The time information will be red colored and is in standby mode.

**Note,** By grabbing the first image automatically the timer will start. When prefer to start the timer manually press the "start timer" button, that's appear when you enable the timer.

**Note,** The images will be shown in your capture window, below the capture window there is a thumbnail overview of the images.

**Note,** The thumbnails size can be modified from:

File→options→General→Thumbnail size.

You also can hold the control key and scroll with your mouse wheel.



#### NOTE

Procedures can also be start from the Procedures menu and for some instruments, from the instrument.



#### CAUTION

When changing procedure, by means of the instrument, take into account that the instrument can be faster than IMAGEnet i-base. Therefore delay capturing until IMAGEnet i-base has the correct procedure in place. You can check the procedure name in the capture session tab.



#### CAUTION


When capturing data from an instrument, the user must not change any instrument or camera settings until the data is downloaded into i-base. Otherwise data properties may not be correct.


### ADD A NEW PATIENT OR SELECT A PATIENT

When starting a new capture session, we need to add patient information to this capture. IMAGEnet i-base patient management window will pop up automatically every time when you start a new capture session.

We can add a new patient or select an existing patient from the patient list.

## ADD A NEW PATIENT

From the patient management window we select the new patient button  from the action bar.

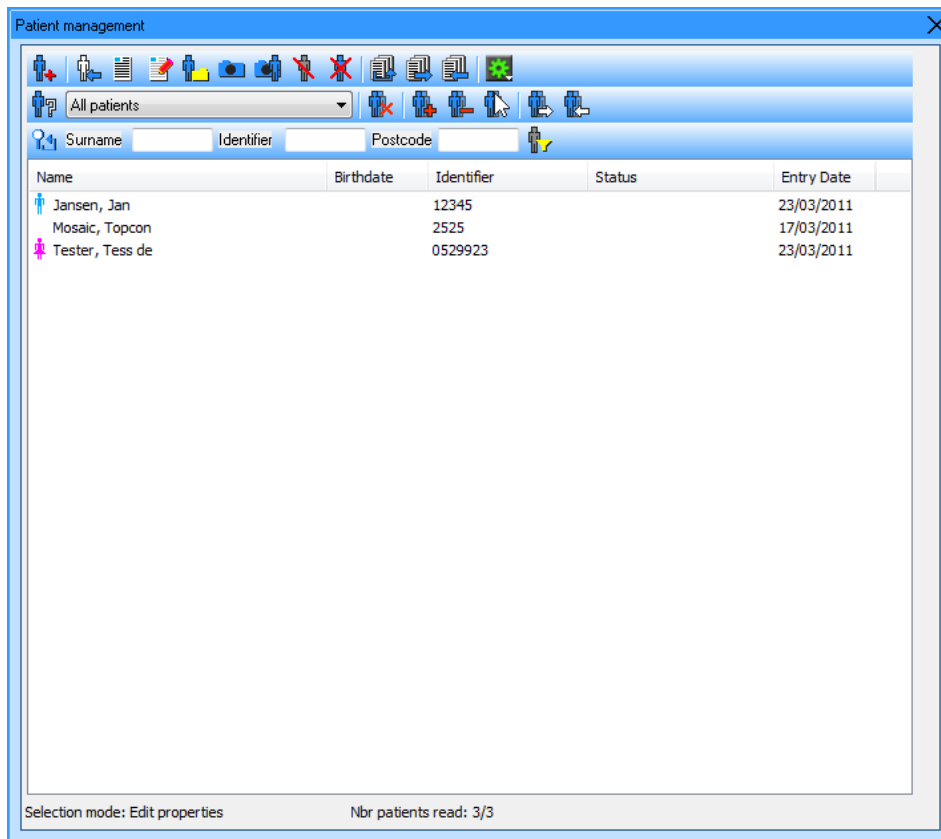


The screenshot shows the 'Patient management' window. At the top, there is a toolbar with various icons, including a blue plus sign icon for adding a new patient. Below the toolbar is a search bar with fields for Surname, Identifier, and Postcode. A table lists existing patients:

Name	Birthdate	Identifier	Status	Entry Date
Jansen, Jan		12345		23/03/2011
Mosaic, Topcon		2525		17/03/2011
Tester, Tess de		0529923		23/03/2011

At the bottom, it indicates 'Selection mode: Edit properties' and 'Nbr patients read: 3/3'.

Fill in all necessary patient information, basic information as Patient ID and Surname are required.

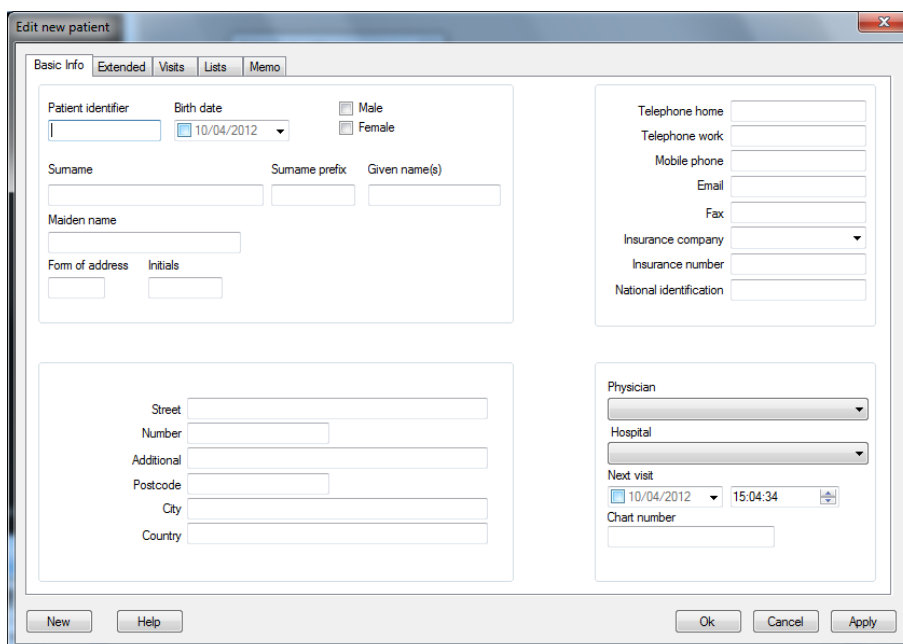


The screenshot shows the 'Edit new patient' window with the 'Basic Info' tab selected. The form contains the following fields:

- Patient identifier:** A text input field.
- Birth date:** A date picker set to 10/04/2012.
- Gender:** Radio buttons for Male and Female.
- Surname:** A text input field.
- Surname prefix:** A text input field.
- Given name(s):** A text input field.
- Maiden name:** A text input field.
- Form of address:** A dropdown menu.
- Initials:** A text input field.
- Telephone home:** A text input field.
- Telephone work:** A text input field.
- Mobile phone:** A text input field.
- Email:** A text input field.
- Fax:** A text input field.
- Insurance company:** A dropdown menu.
- Insurance number:** A text input field.
- National identification:** A text input field.
- Street:** A text input field.
- Number:** A text input field.
- Additional:** A text input field.
- Postcode:** A text input field.
- City:** A text input field.
- Country:** A text input field.
- Physician:** A dropdown menu.
- Hospital:** A dropdown menu.
- Next visit:** A date and time picker set to 10/04/2012 15:04:34.
- Chart number:** A text input field.

At the bottom, there are buttons for 'New', 'Help', 'Ok', 'Cancel', and 'Apply'.

Press OK to return to the capture session window.

Patient information is now also available in your capture session window.

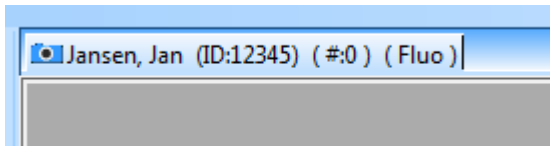

Patient name = Jan Jansen


Patient ID = 12345

"#" amount of captured images

Fluo = procedure name

## SELECT A PATIENT

The patient management window will start automatically when a new capture session starts. Now we will select a patient from the patient list. From the query bar we can filter patient on surname, patient ID and postcode. To start the capture for this patient double click on the patient name.

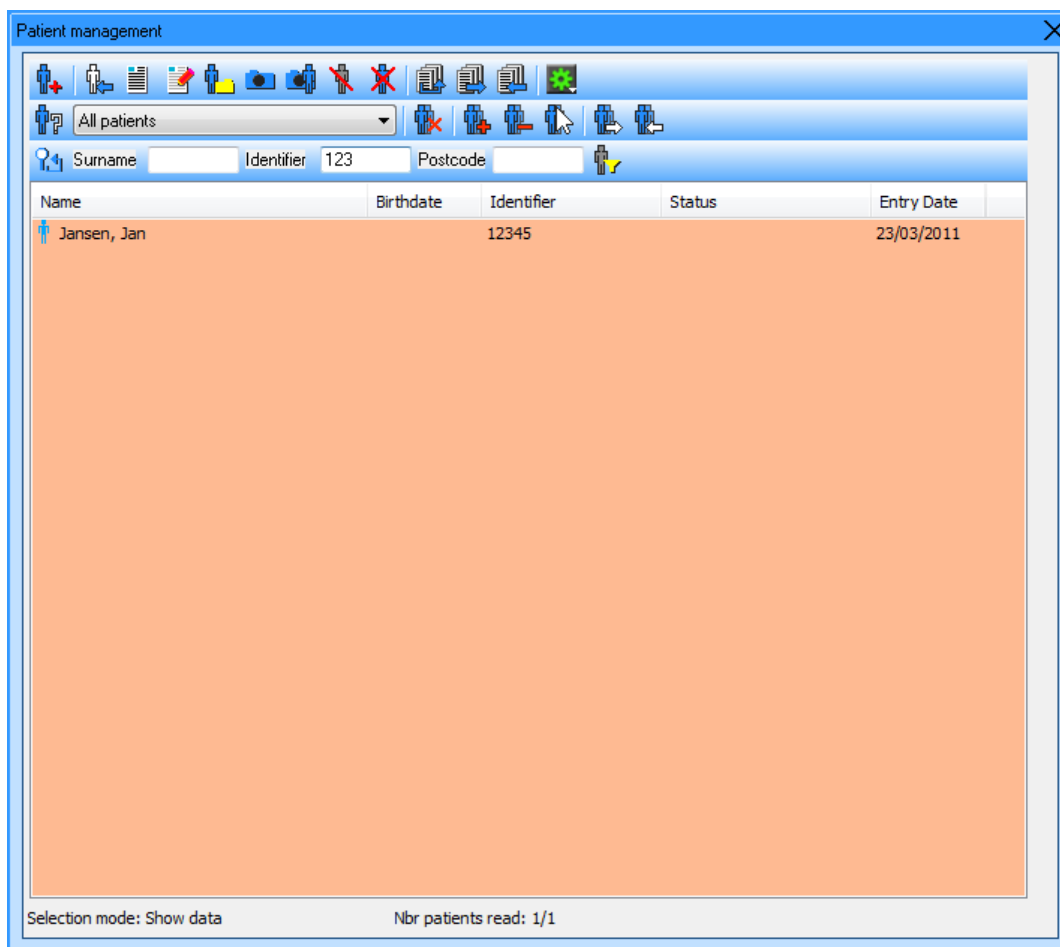


### NOTE

The background of the patient list will be colored Tan when a query is active.



## CAPTURE WINDOW READY TO CAPTURE IMAGES

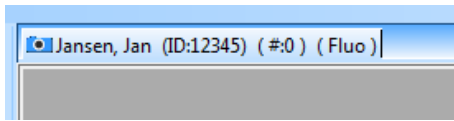


Check if procedure and patient information are correct when capture session window is open.  
Now you are ready to start grabbing images.


## START TIMER

To enable the timer press the “Enable timer” button.  
The timer information will be red colored and is in standby mode.

By grabbing the first image automatically the timer will start.  
When you prefer to start the timer manually press the “Start timer” button.  
“Start timer” button appear when you enable the timer.

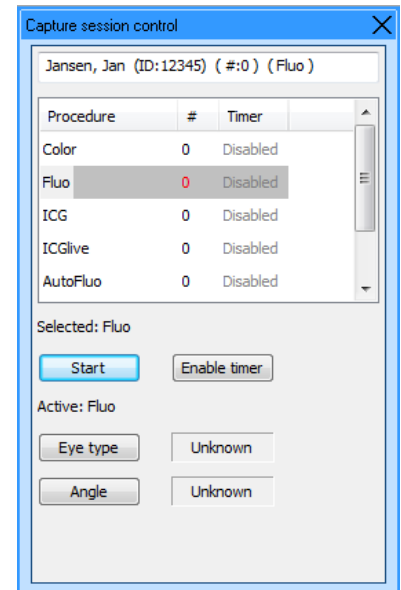


## DELETE IMAGES DURING CAPTURE

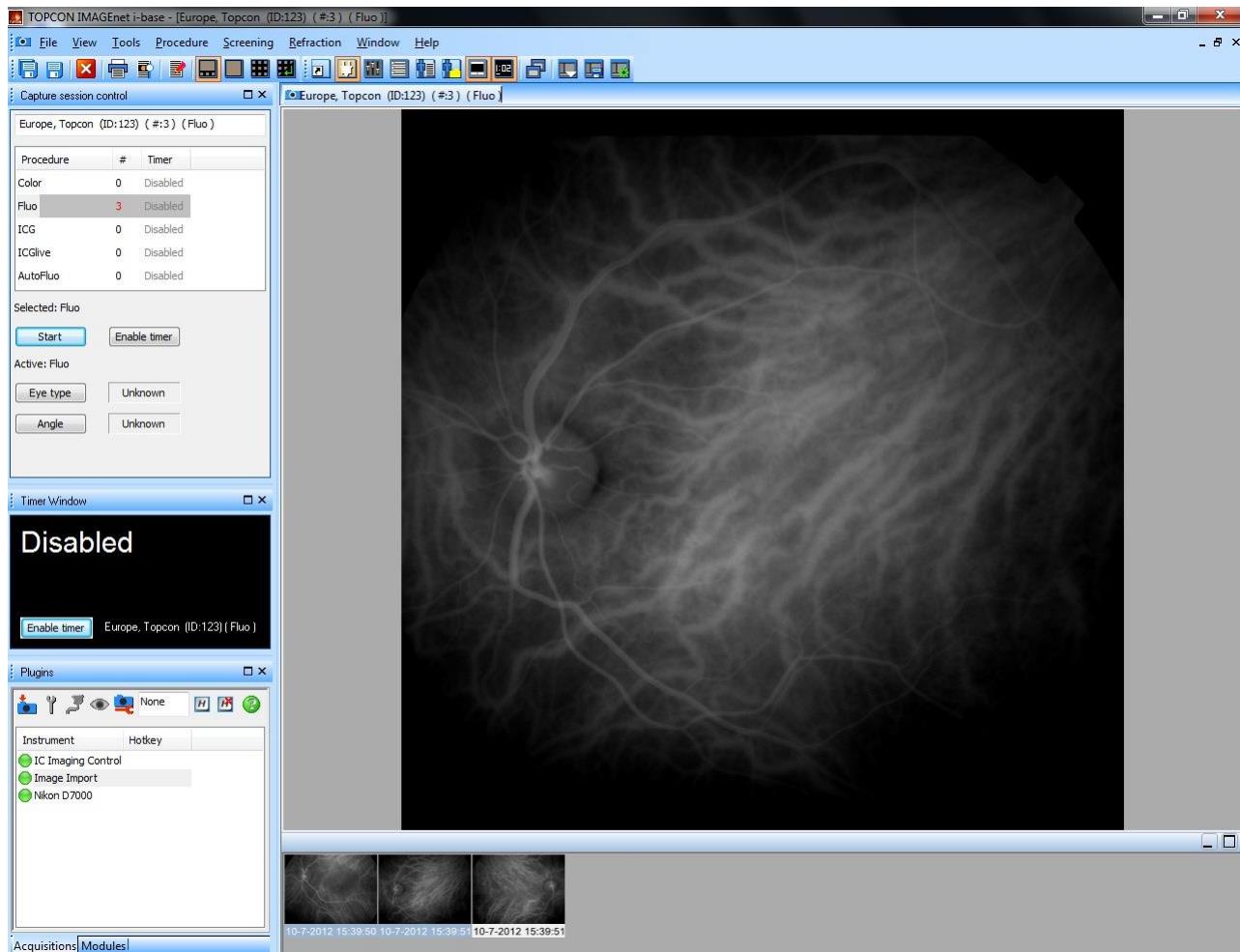


During the capture you can scroll back through the thumbnails and view the captured images by making the thumbnail active. Images that are not correct you can delete by selecting 1 or more (with shift or ctrl) thumbnails and use the delete key or with the right mouse on a thumbnail and select "delete selected"

## UNDELETE IMAGES DURING CAPTURE

Restore of deleted images can be done through the shortcuts window → acquisition and use the recycle bin button "show/hide deleted".

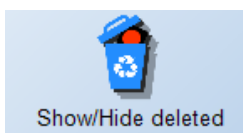


Right mouse on the thumbnail and use undelete option.

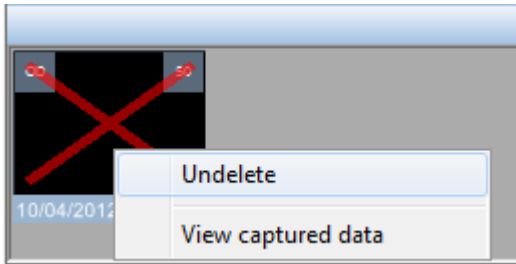
**NOTE**

Images that are deleted can only be restored when the capture session is active

## SAVE CAPTURED IMAGES

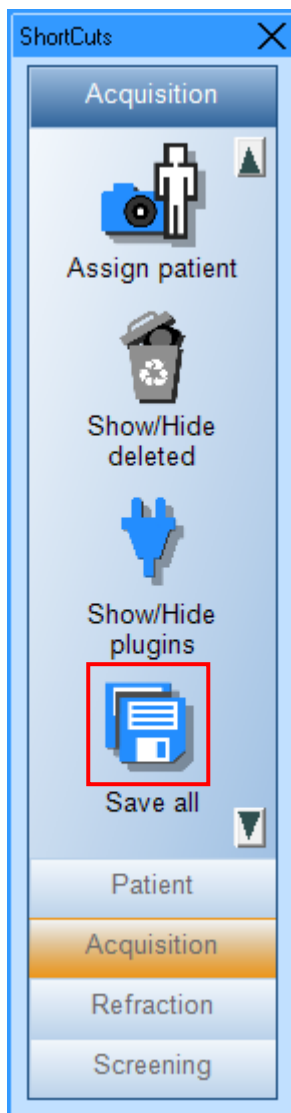

To save all the captured images you can use the “Save all” button from the shortcuts window or the “Save all” button in the standard tool bar

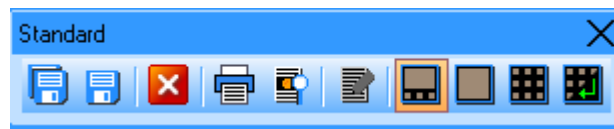

Button to save all images




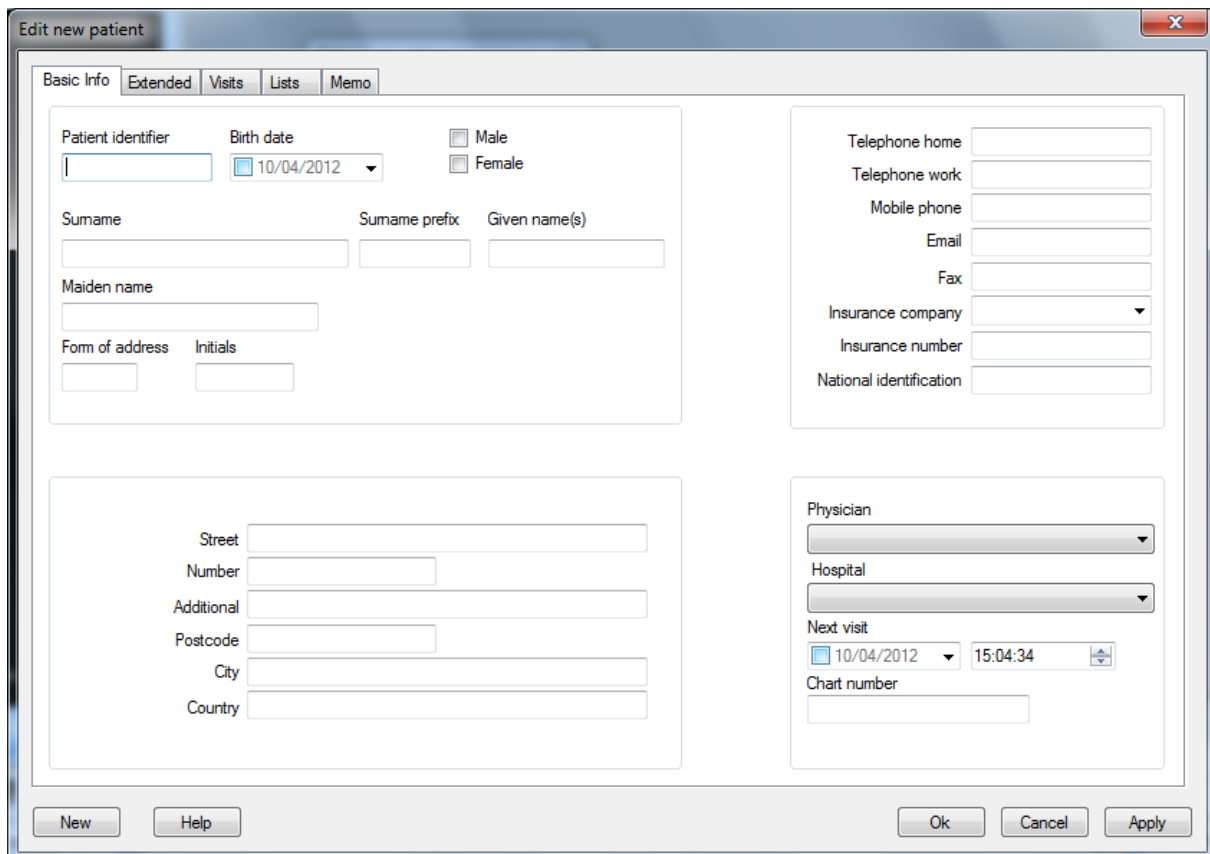
Button to save selected images

## PATIENT MANAGEMENT / DATA

From the shortcuts bar → Patient, you can add a new patient, edit the patient and show the patient data.

### ADD NEW PATIENT

From the patient management window we select the new patient button  from the action bar or from the file menu from the windows tool bar select “new patient” and fill in all necessary patient information.



Basic information as Patient ID and surname are required, press “OK” and the patient will be added to the patient list and the patient management window will close.

Press “New”: Your patient is added to your patient list, but you are able to continue adding new patient.


Press “Apply”: Your patient is added to your patient list, now you can add also other information from the different tabs in the patient properties window.



#### NOTE

There is also a minimal patient properties window available.

To use this, select in the option menu → general → patient management “minimal patient dialog”.

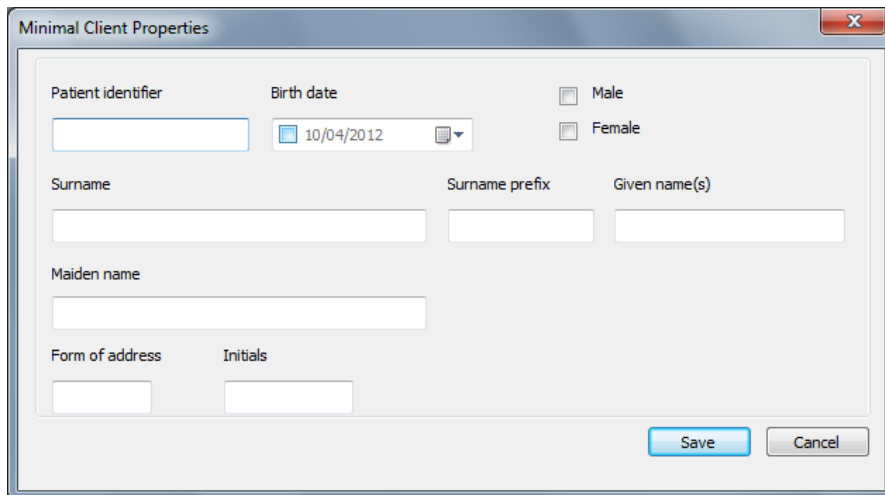

**Minimal Client Properties**

Patient identifier:  Birth date:   ☐ Male ☐ Female


Surname:  Surname prefix:  Given name(s):


Maiden name:

Form of address:  Initials:

## EDIT (MANAGE) PATIENT INFORMATION

From the shortcuts bar select “Manage Patient” , now you are in the edit patient mode.

Select a patient and press the “edit properties” button  or double click on the patient name. The edit properties window will appear, now you are able to edit patient information. Press the OK button to save your changes.

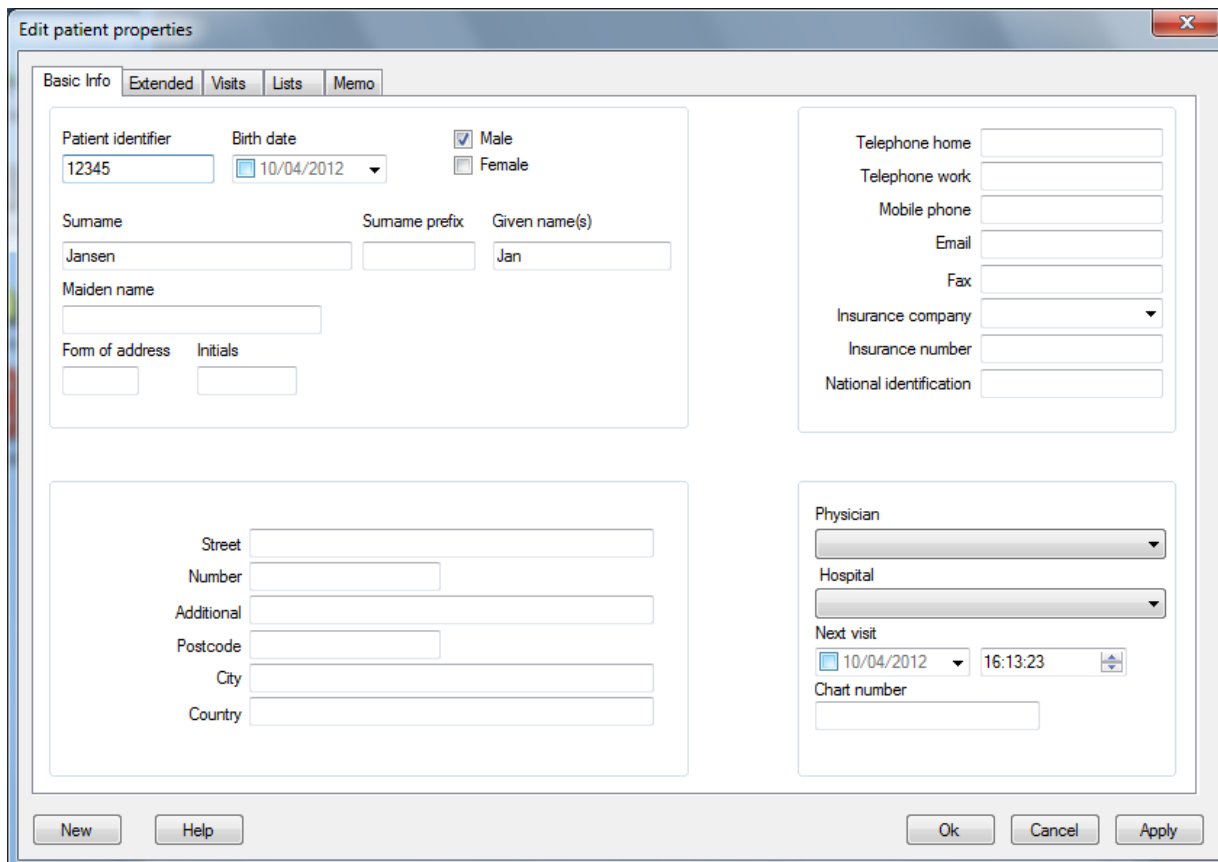

**Edit patient properties**

Basic Info | Extended | Visits | Lists | Memo

Patient identifier:  Birth date:   ☒ Male ☐ Female

Surname:  Surname prefix:  Given name(s):

Maiden name:

Form of address:  Initials:


Telephone home:   
Telephone work:   
Mobile phone:   
Email:   
Fax:   
Insurance company:   
Insurance number:   
National identification:


Street:   
Number:   
Additional:   
Postcode:   
City:   
Country:

Physician:   
Hospital:   
Next visit:    
Chart number:

## SHOW PATIENT DATA

From the shortcuts bar select “show patient data”


Select a patient and press the “show data” button  or double click on the patient name  
The patient data window will shown, now you are able to view patient data.

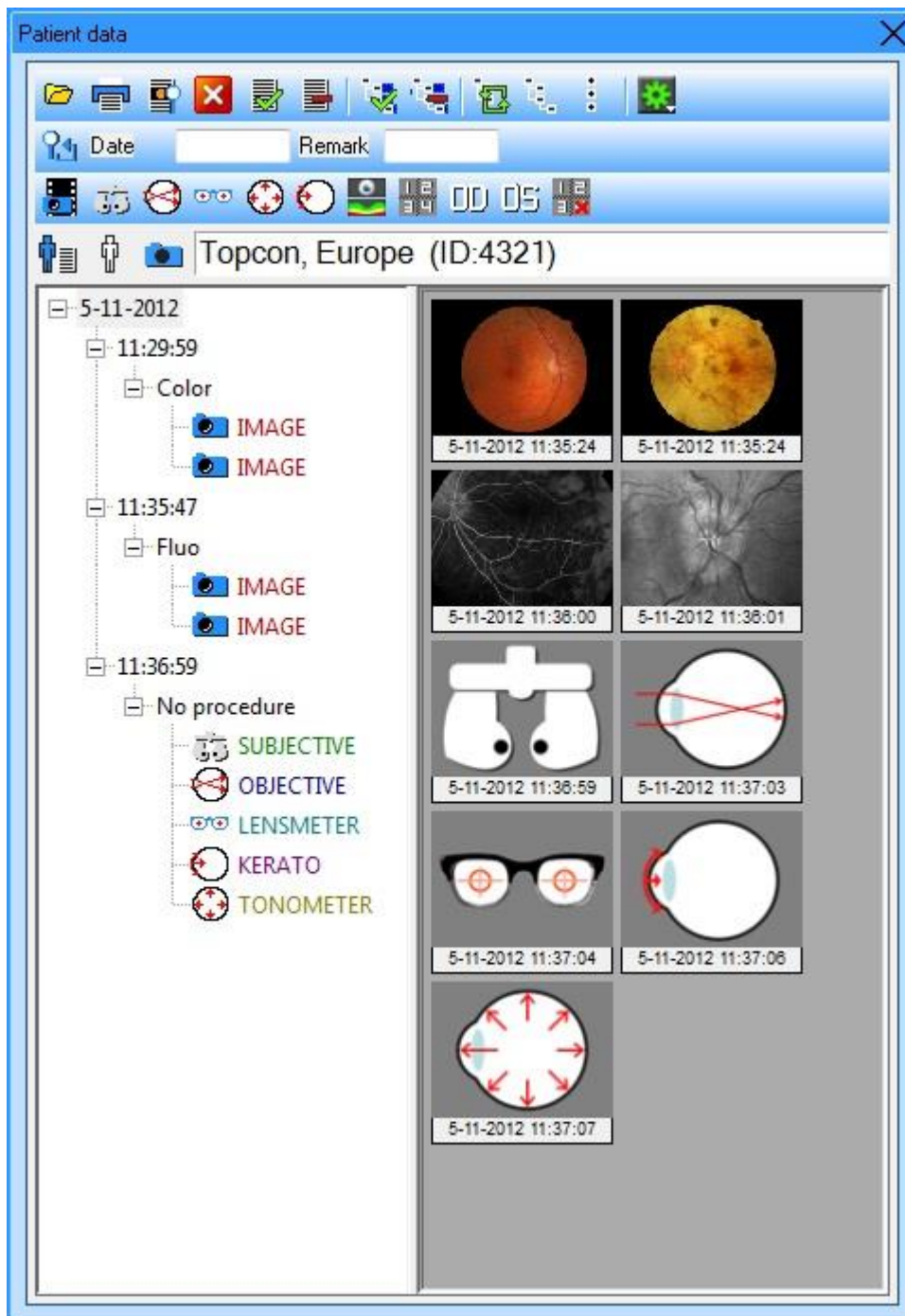



### NOTE

From the patient management window you can also select new, edit, show data etc. with the right mouse button.

## OPEN PATIENT DATA


To open data or images use the open data button  or select a visit or procedure and use right mouse button to select open data. If you only want to open one image, you can double click the image.

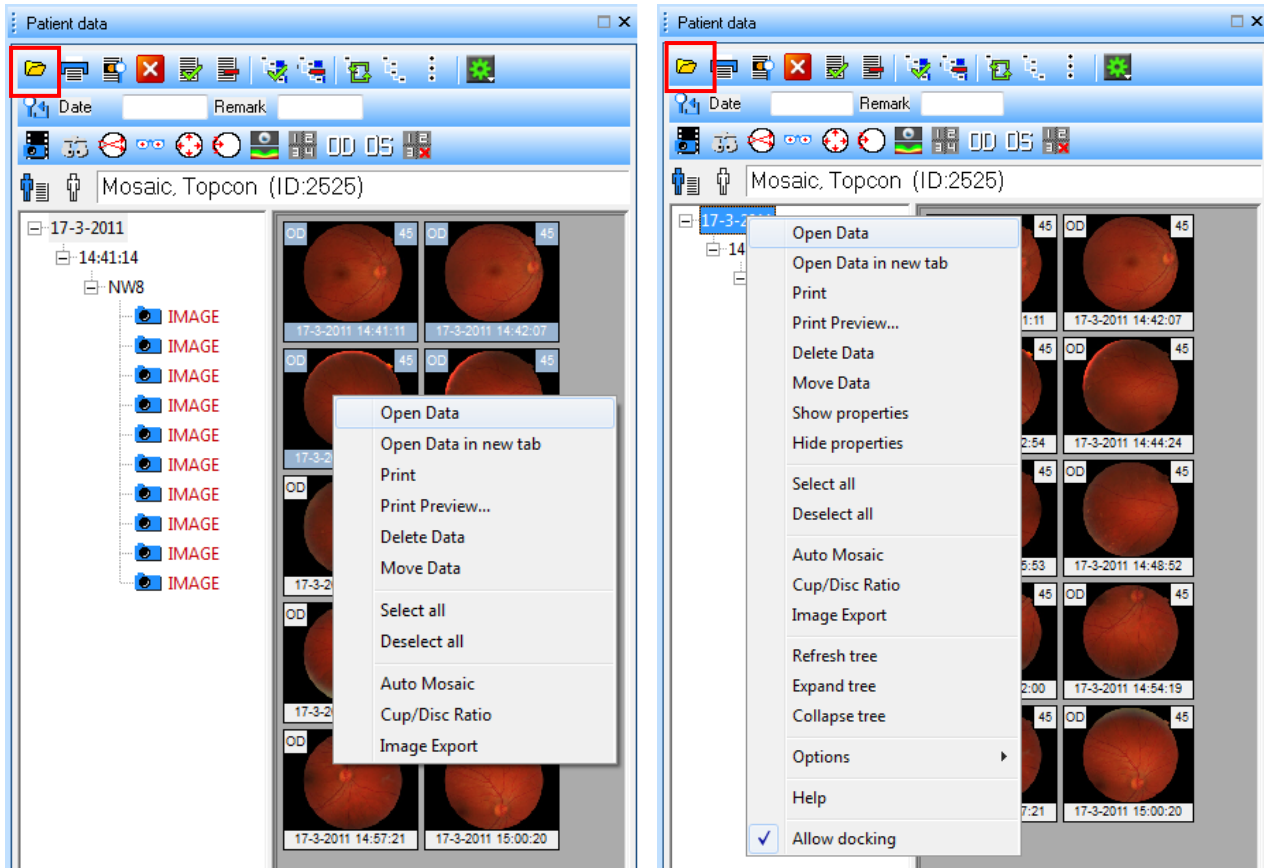


You can also make a selection of images (ctrl + left mouse button) and use the open data button  or use right mouse button to select open data.

### NOTE




To close the patient data window after opening images or data, select in options  from the patient data window "Close window after opening items" or use the right mouse button and go to options.

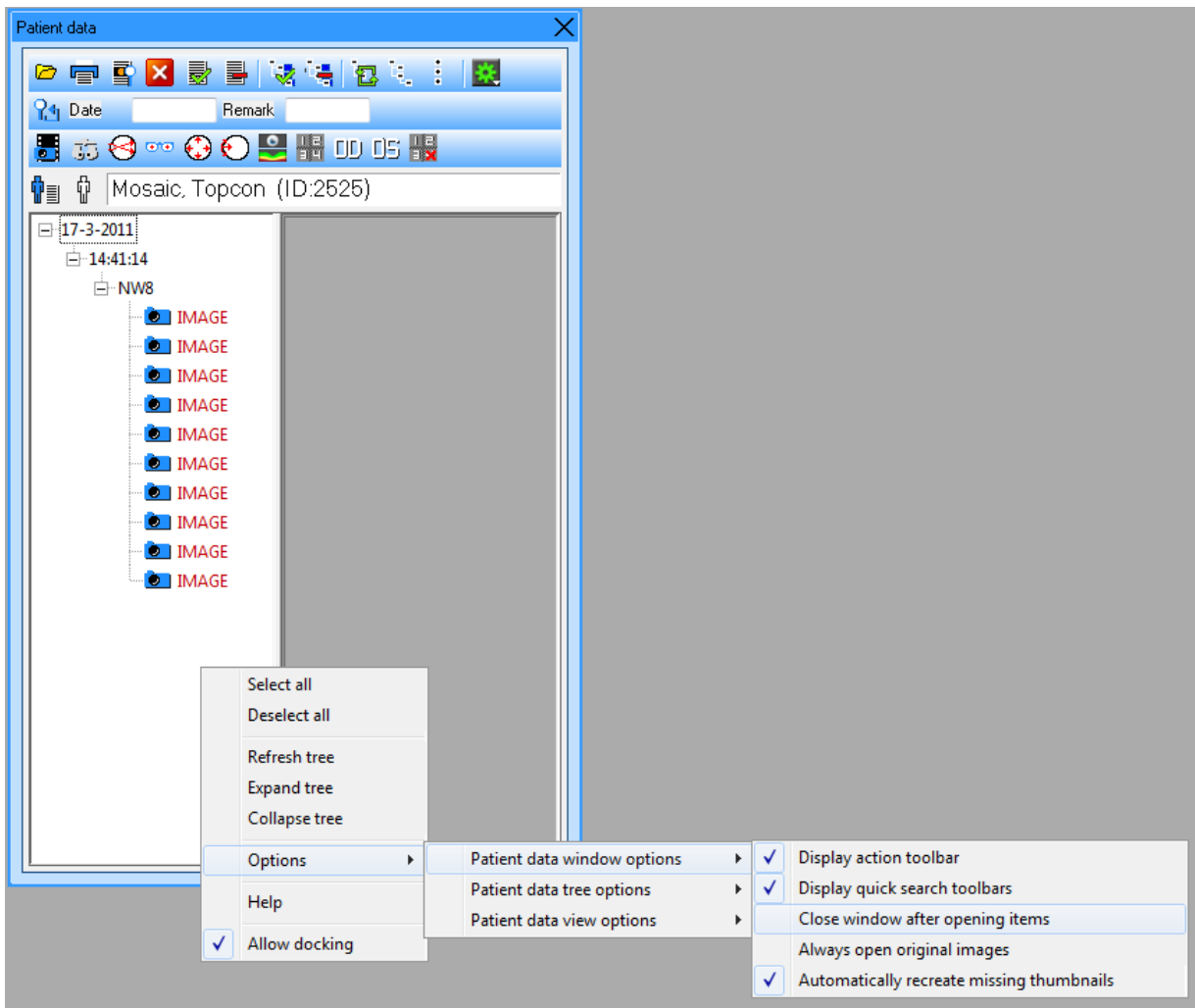

### NOTE

✓ To open data or images directly in a multi-up window, select from File > options > general > opened data, "Present opened data always in a multi-up".

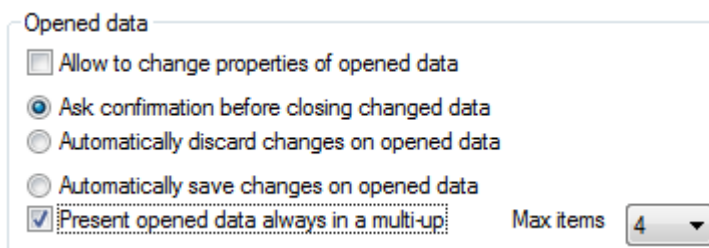


## WORKING IN A MULTI-UP

Open data in a multi-up view is by default disabled, you can work manually by opening images and double click on the first image, a 4 up will be displayed. By selecting images from the thumbnail bar you can add, removed or swap (hold ctrl key) images.


By double click in the 4-up you will be switch again to your single-up



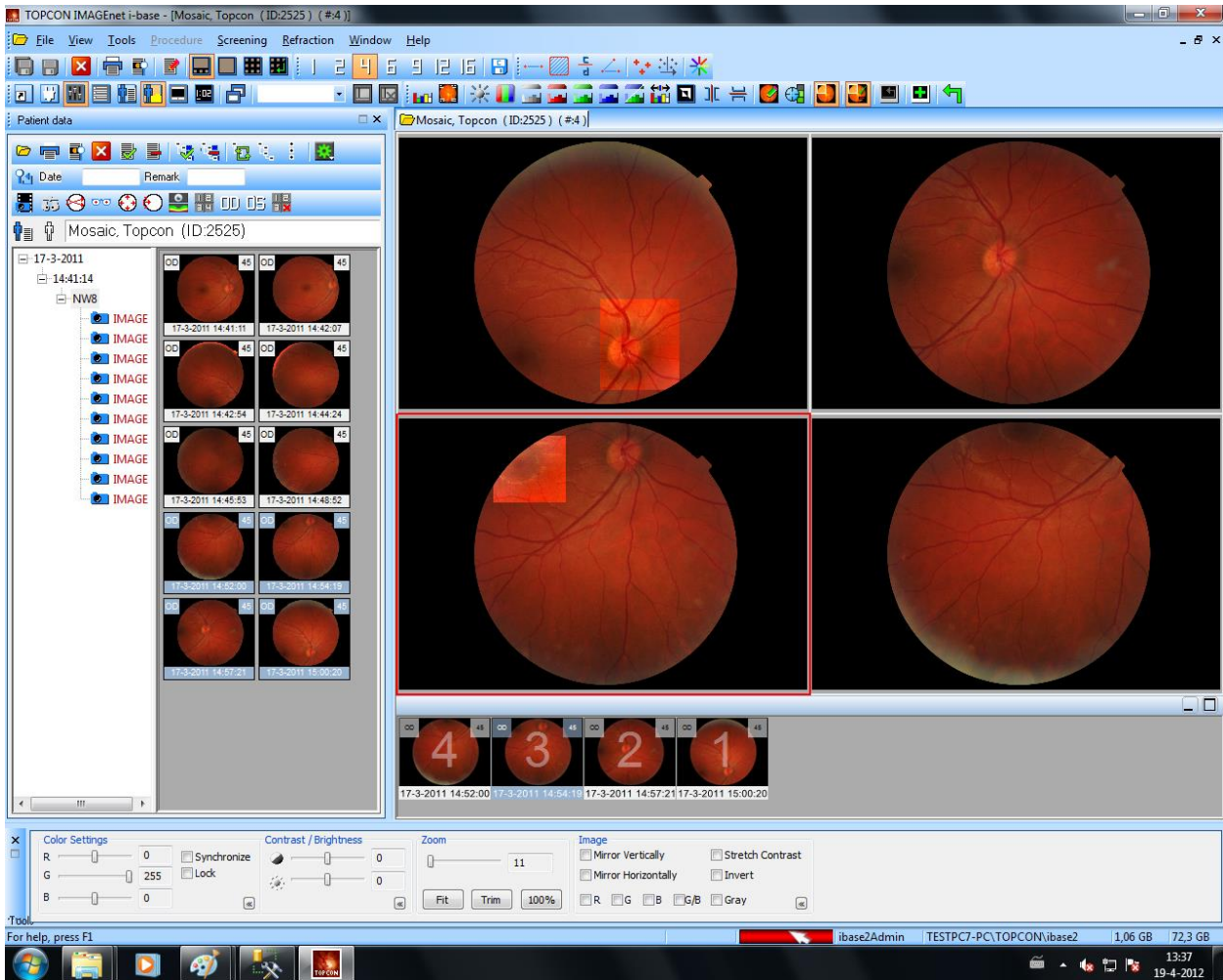


## COMPOSITION TOOL BAR

The composition toolbar provides fast access to several composition functions.

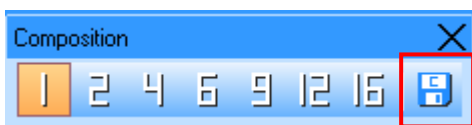

To save the current views in the selected acquisition or viewing tab as a composition select the "Save button". A dialog pops up where you can enter the composition name, press OK to Save

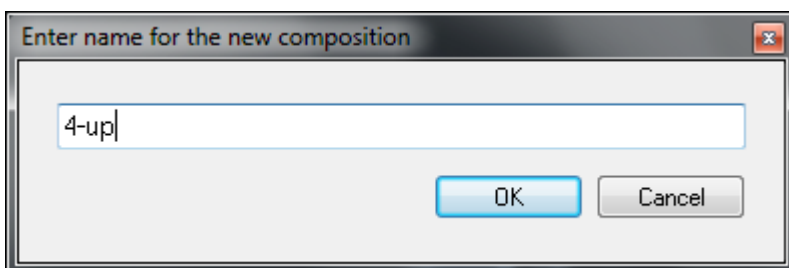







**NOTE**

To open data or images directly in a multi-up window, select from File > options > general > opened data, "Present opened data always in a multi-up".


## FITTING IMAGES BY ZOOM

The image zoom toolbar provides fast access to several zooming functions.



	Normal size:	zoom the image or movie to the original size.
	Trim zoom:	zoom the image so that the black background is hidden.
	Zoom fit:	zoom the image or movie to fit the view.
	Rectangle zoom:	when this option is enabled, you can draw a rectangle in the image to zoom in to.
	Zoom glass:	with this option you get a zooming glass, which you can use to zoom into the image (this option is only available in the EyeDoctor and EyeViewer modules).
	View pan window:	shows a window with a thumbnail of the currently viewed image and the current 'Zoom rectangle'. Within this 'Pan window', you can drag the 'Zoom rectangle' with the mouse to pan over the image.

To fit open your images as default in your window fully visible you have to set some media options.  
Select File → Options → Media .

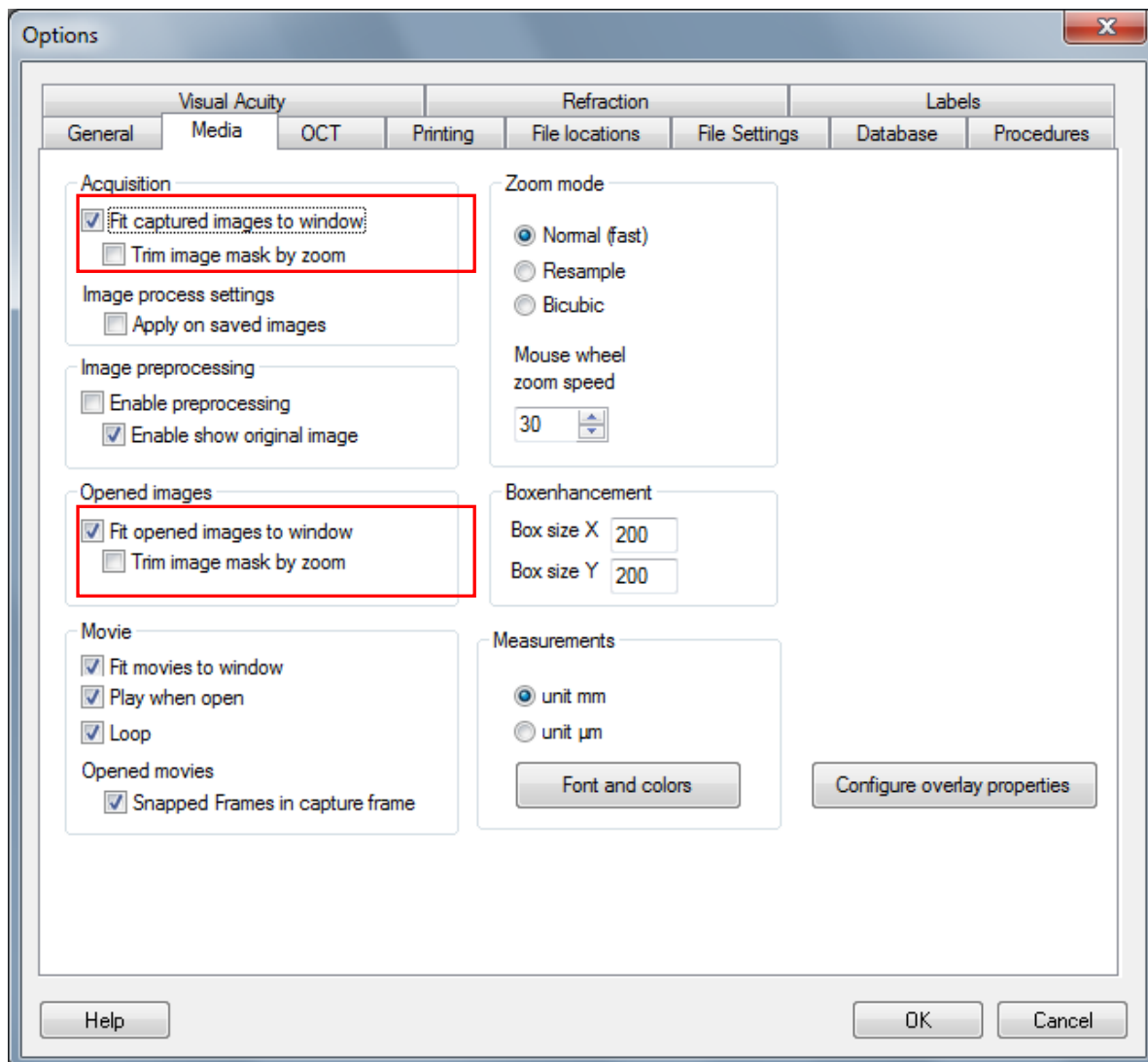

## IMAGE PROCESS

When select a image you can adjust manually or automatically with help of the Image process bar.



Shows color histogram



Shows Image properties (you also can use the F8 key)



Contrast and brightness adjust



Adjust color balance



Separate colors- 8 bit grey, red, green and blue free.



Shows red free image



Auto stretch contrast



Invert



Mirroring



Annotate



Sharpening



Box enhancement



Show / hide box enhancement



Open original image




Apply adjustments in new image



Undo

## TOOLS WINDOW

To show this tool window please go to view and then select tools.

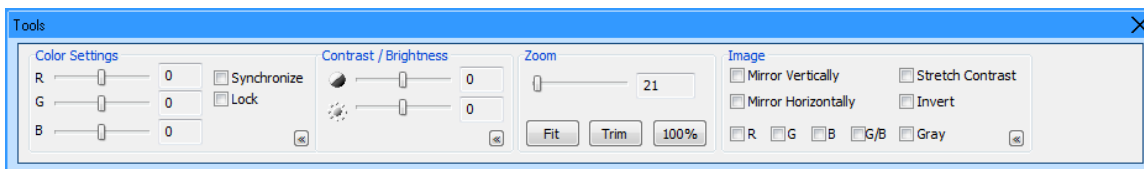

This tool can be used as an easy and fast way to process opened or captured images. For captured images this tool can only be used for viewing. Captured images can never be changed. For opened images this tool can be used for viewing, but if needed you can also apply the changes in a new image with 'Apply in new image' from the Tools menu or with the 'Apply in new image' button in the Image process toolbar. Now a new image is created and saved. Original captured images can never be changed.



## PLUGIN WINDOW

From the plugins window the installed plugins can be managed.

Two types of plugins are available within IMAGEnet i-base: 'Acquisition' plugins and 'Module' plugins. The two types are separated in the plugins window with two tabs. When you select a tab you get a list of the installed plugins of that type.

Acquisition plugins can send data from external sources (for example a digital camera) to IMAGEnet i-base.

Module plugins can receive data and patient information from IMAGEnet i-base and optionally send data back to IMAGEnet i-base.

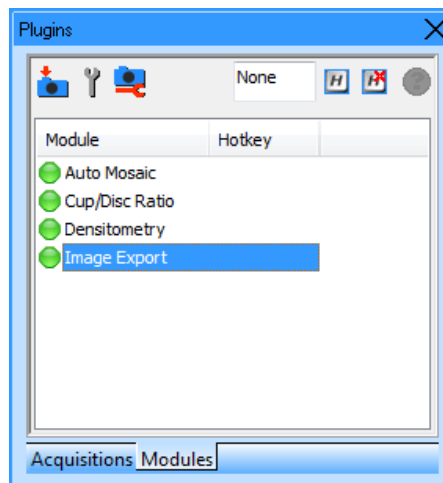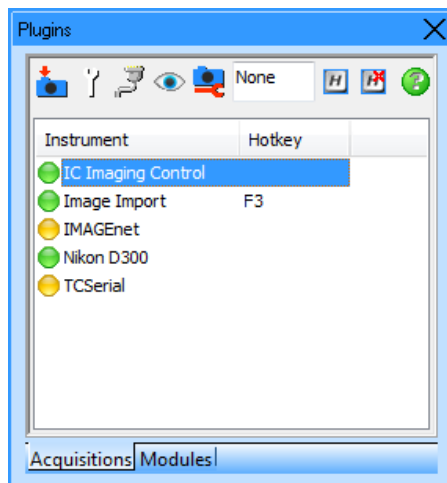

**Note,** A plugin can have three statuses:


green: available  
yellow: not available  
red: not loaded

start the selected plugin; for acquisition plugins this usually means trigger the camera; for module plugins this usually means start the module with the selected data items.


To setup a plugin just right-click on it.

Modules can start by:


*Double click in plugin window*

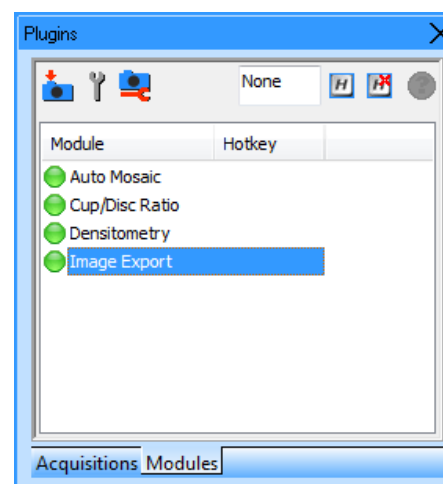

*Selecting the thumbnail(s) in patient data window and right click*

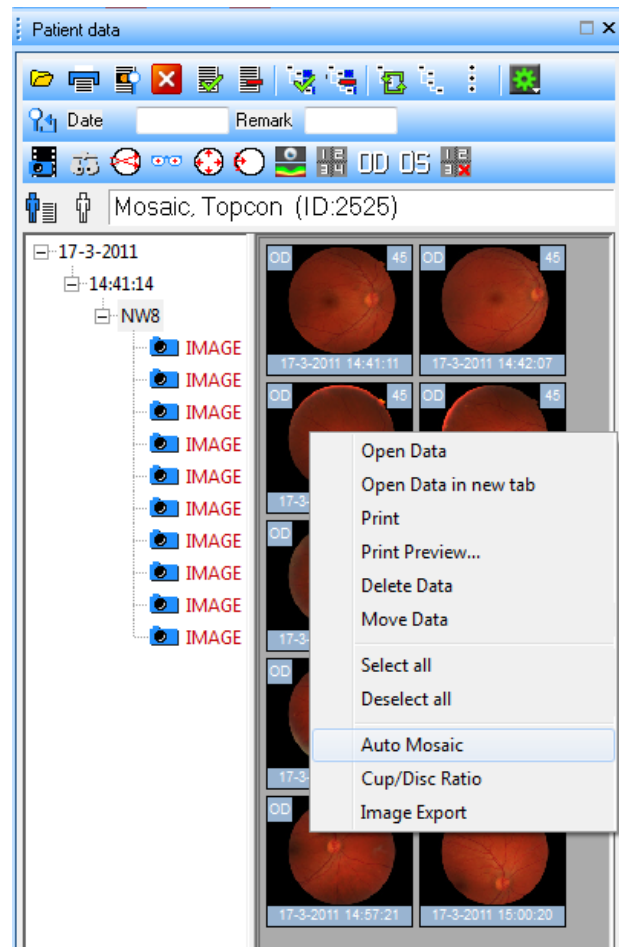


*Selecting the thumbnail(s) of open image(s) and right click*

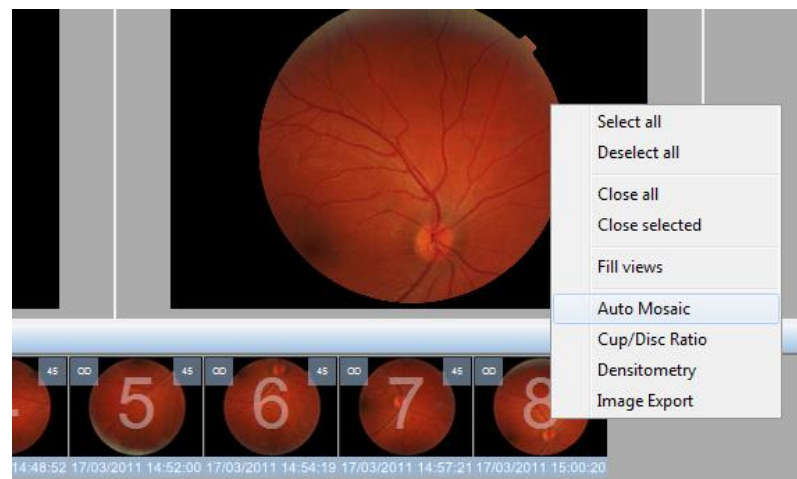


## MODULES

### EXAMPLE AUTOMOSAIC

Select all your mosaic images in the Patient data window and right click on one image, now select Automosaic and the Automosaic window will open.

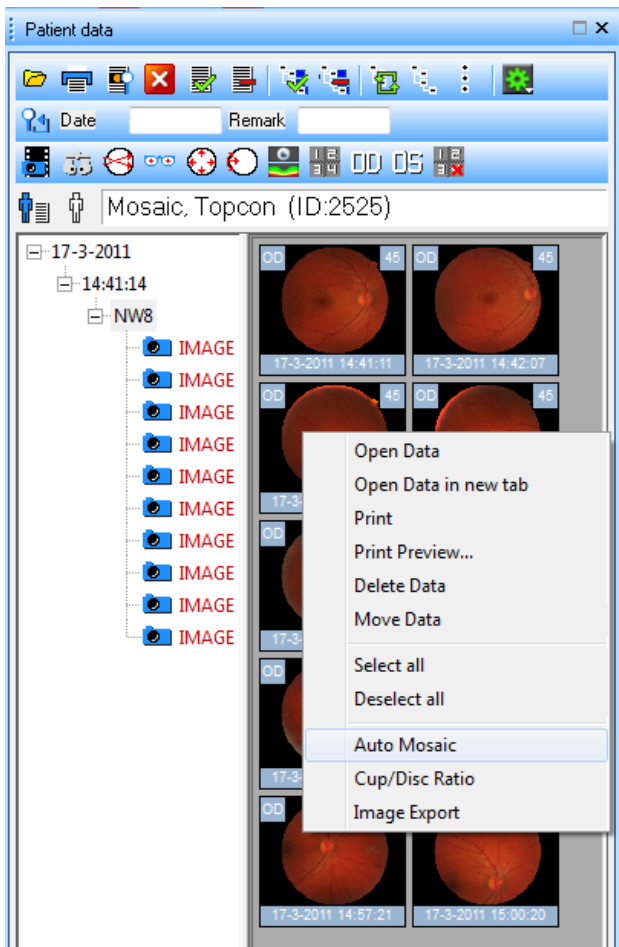


If there is different information about eye and angle in the selected images, this screen will open first:

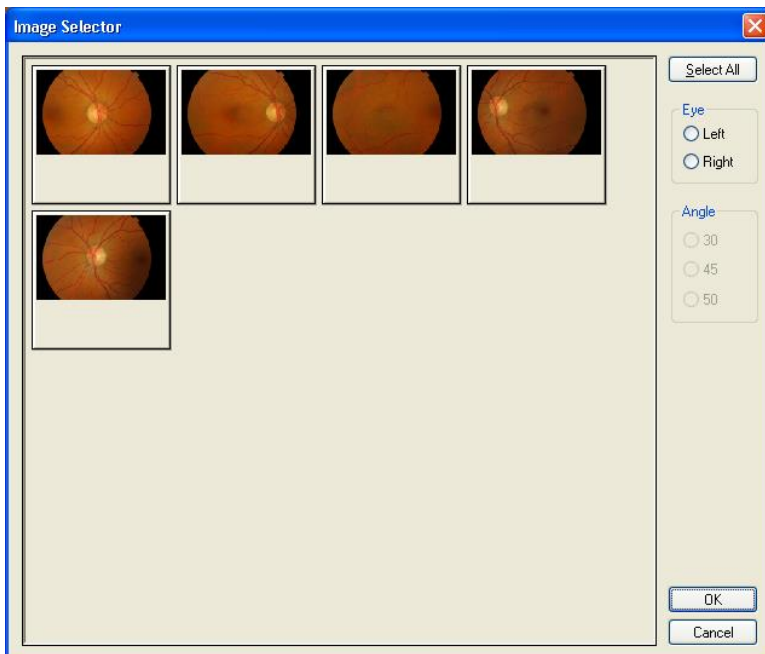

Here you can select eye and angle. Only the images that meet the selection are going to the Automosaic module.

## EXAMPLE IMAGE COMPARE

If you want to use Image Compare, you first have to make some registration points on the 2 images. You have to make 3 registration points at least and save the images with the registration points.

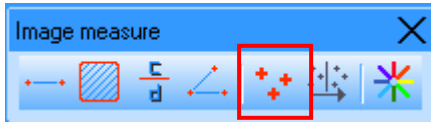


registration points you will find on the Image measure toolbar.




After saving the images you have to select the 2 thumbnails and right-click on them and select image compare or double-click image compare in the plugin window.



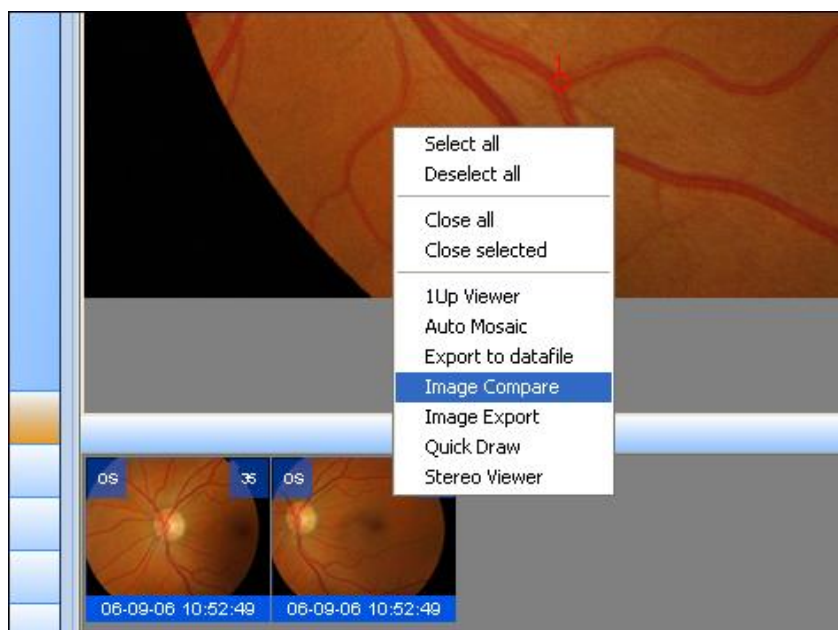


Now this window will open and you can start comparing:

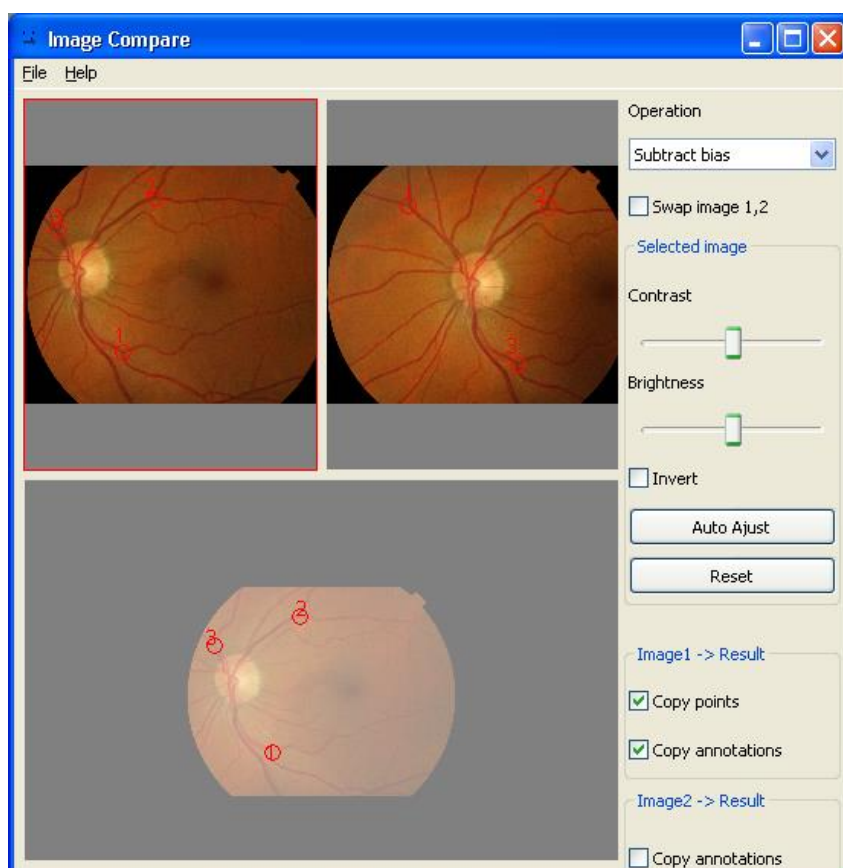

## EXAMPLE DATAFILE


To Import a DataFile please double-click the DataFile module in the plugins window or right-click and select start. Now this window will open:

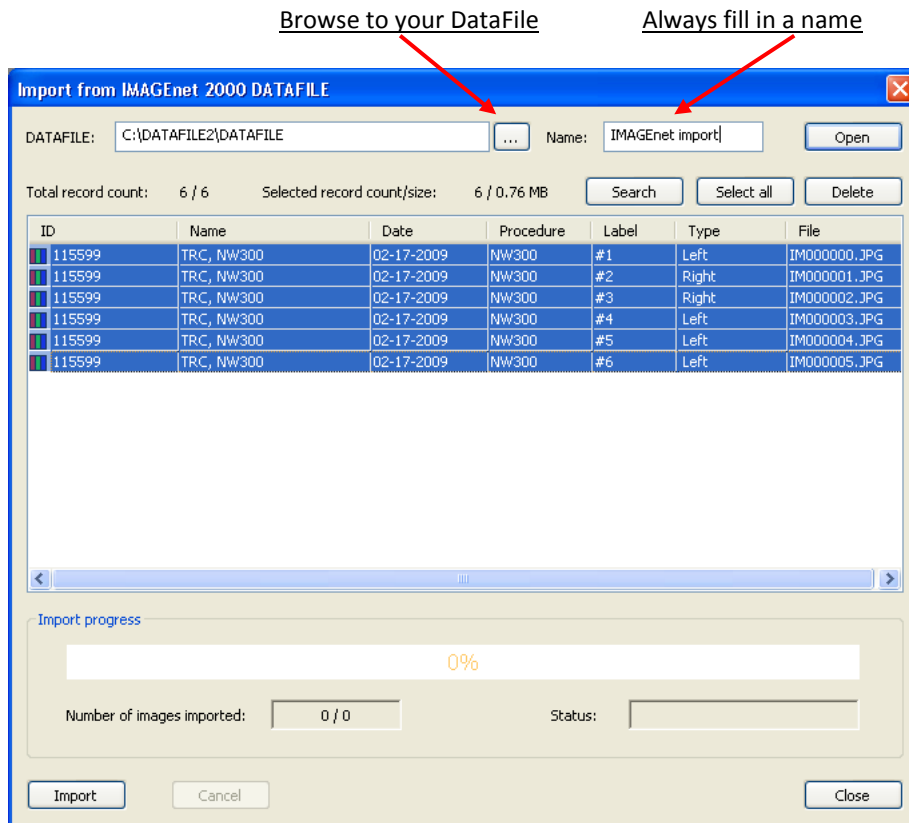

Now click Import and your Selected Patients and images will be imported.

## MEASUREMENTS

To make a measurement on an Image you select the measurement toolbar:

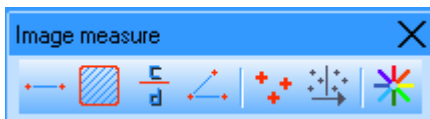

Line Measurement



Area Measurement



C/D Ratio Measurement



I-Torsion Measurement

**NOTE**


The Image measure toolbar is only active with opened data.

## PRINT PATIENT DATA

The print function to print patient data and/or properties within IMAGEnet i-base is based on templates. A default portrait and landscape template are installed.

There are several ways to print patient data:


*From the File menu*

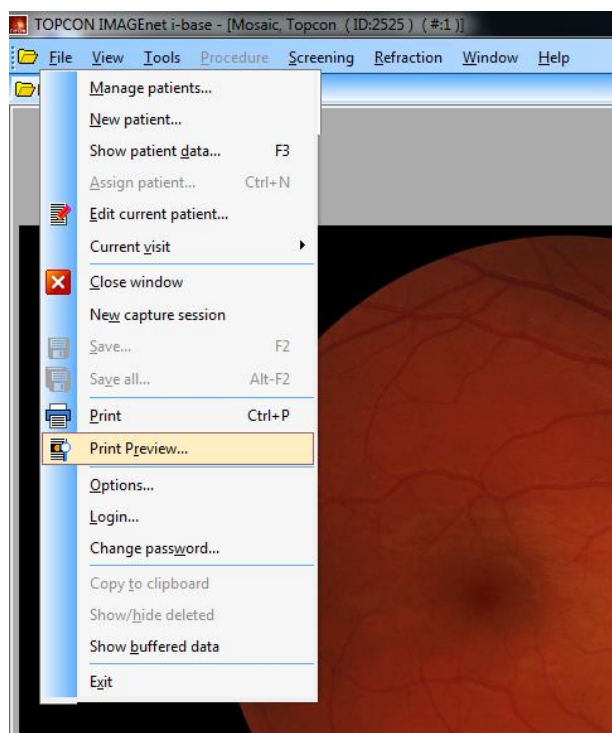

*From the standard toolbar*



From the Patient data window

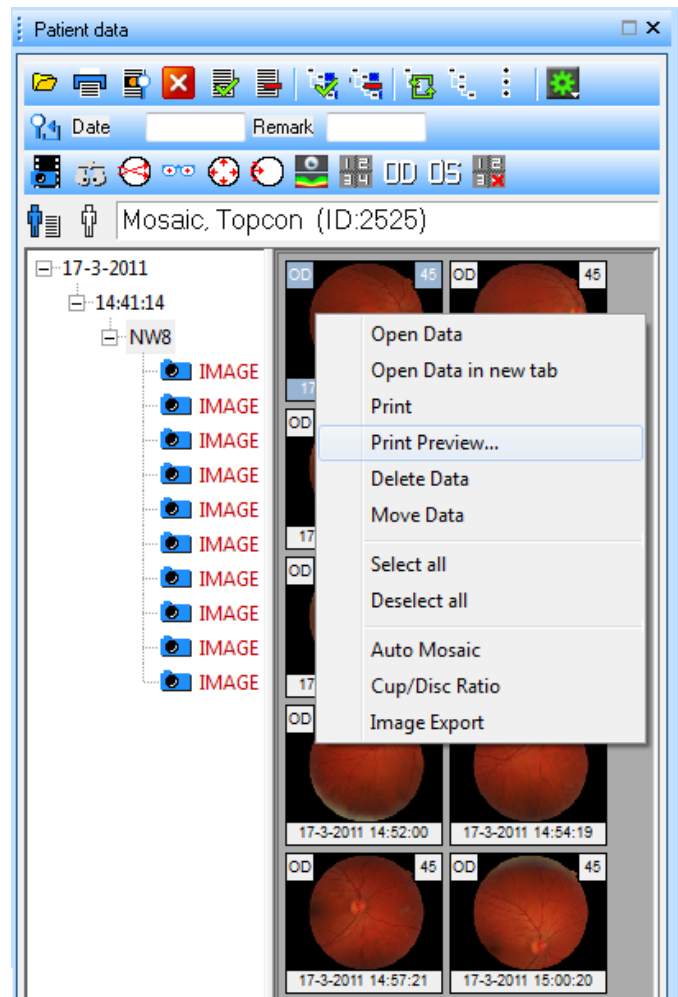

## PRINT PREVIEW

Print Preview shows a print preview of the selected acquisition or viewing tab.

The preview window is setup with the default print settings as set in printing options. These settings can be changed in this window.



Mask image: Mask the selected images



Contrast stretch image: Stretch the selected images




Trim zoom: Trim the selected images




Show/hide properties: Put the overlay properties in the selected images



Configure "show/hide properties": Set here the overlay properties. Use the button 'Overlay' to put the overlay properties in the selected images.

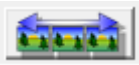

Select all the images in the thumbnail list

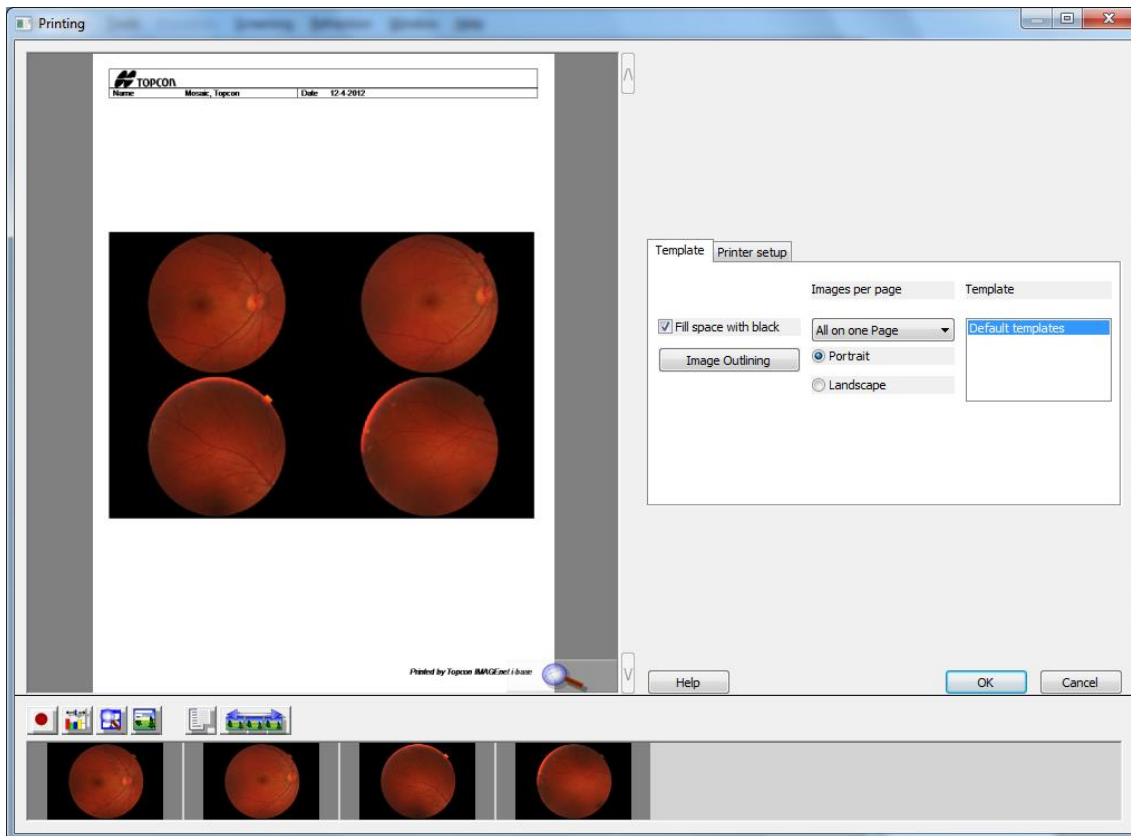

Press OK to print, print preview will closed after print task is given.



TOPCON EUROPE MEDICAL BV  
11 - 2908 LJ CAPELLE AAN DEN IJSSEL  
T: +31 (0)10 4585077  
E: [medical@topcon.eu](mailto:medical@topcon.eu)  
Website: [www.topcon-medical.eu](http://www.topcon-medical.eu)